



# EVALUATION OF CIGUATOXIN-INDUCED CGRP RELEASE IN MOUSE AND RAT SKIN PREPARATIONS

## QUANTIFY CGRP LEVELS IN RODENT SKIN PREPARATIONS USING THE BERTIN BIOREAGENT CGRP ELISA KITS

Filip Touska, Simon Sattler, Philipp Malsch, Richard J Lewis, Peter W Reeh, Katharina Zimmermann  
*Department of Anesthesiology, Friedrich-Alexander-University Erlangen-Nuremberg, University Hospital Erlangen*

### CONTEXT

Ciguatoxins are a group of toxins that are present in certain fish and are responsible for a type of severe food poisoning called ciguatera. The consumption of tropical and subtropical fin fishes contaminated by ciguatoxins can cause several painful neurological, gastrointestinal, and cardiovascular symptoms. Ciguatoxins are known to be powerful voltage-gated sodium channel activator toxins, which partially explains some of the sensory symptoms associated with ciguatera. However, for many of the neurological symptoms caused by ciguatoxins, the underlying biological mechanisms remain poorly understood. It has been recently shown that ciguatoxins are extremely effective at releasing calcitonin-gene related peptide (CGRP) from nerve terminals. **(1)**. Here, we use the **CGRP ELISA kit** to evaluate the release of CGRP in mouse and rat skin preparations induced by Pacific Ciguatoxin-1 (P-CTX-1).

### MATERIALS & PROTOCOLS

- **Skin preparations:** The skin flaps of the lower hindpaws of young male Wistar rats (70–80 g; n = 39) and male C57BL/6J mice (15–35 g; n = 136) were used for the release experiments. The preparation of the skin started at knee level and excluded the toes and it spared larger vessels, and saphenous and peroneal nerve stems. The obtained murine skin flaps had an average weight of 0.10 g while the average weight of rat skin flaps was 0.28 g. They were wrapped around acrylic glass rods with the corium side exposed to the surrounding solution and were fixed with surgical silk threads. The mounted skin flaps were then washed for 30 min in carbogen-gassed synthetic interstitial fluid (SIF, pH 7.4) whose composition is described in **(1)**.



# EVALUATION OF CIGUATOXIN-INDUCED CGRP RELEASE IN MOUSE AND RAT SKIN PREPARATIONS

## MATERIALS & PROTOCOLS

- **Stimulation procedures:** After the initial 30 min wash-out in SIF, the skin flaps were placed for 5 min each into pre-warmed glass tubes mounted in a shaking bath at 32 °C. One day prior to the experiments, the tubes were treated with Sigmacote© (Sigma-Aldrich, Taufkirchen, Germany) to prevent adhesion of CGRP and P-CTX-1 to the glass surface. Sigmacote® was refreshed after each 10 experiments. All tubes were filled with 0.7 mL of SIF or stimulation solution containing various concentrations of P-CTX-1 or drugs
- **ELISA assays:** iCGRP (immunoreactive calcitonin gene-related peptide) content in the eluates was measured using the **CGRP rat ELISA kit** (A05482, Bertin Bioreagent, France) as described in (2,3). Briefly, the samples were mixed with 5-fold concentrated commercial CGRP enzyme immunoassay buffer (200 µL sample + 50 µL buffer). The enzyme immunoassays were run on 96 well plates, which were photometrically determined using a microplate reader (Dynatech, Channel Islands, UK). The minimum detection limit observed was 5 pg/mL. The intra and inter-assay coefficients of variation with repeated measurements were 10–15%. Results can be seen in **Figure 1**.

## RESULTS

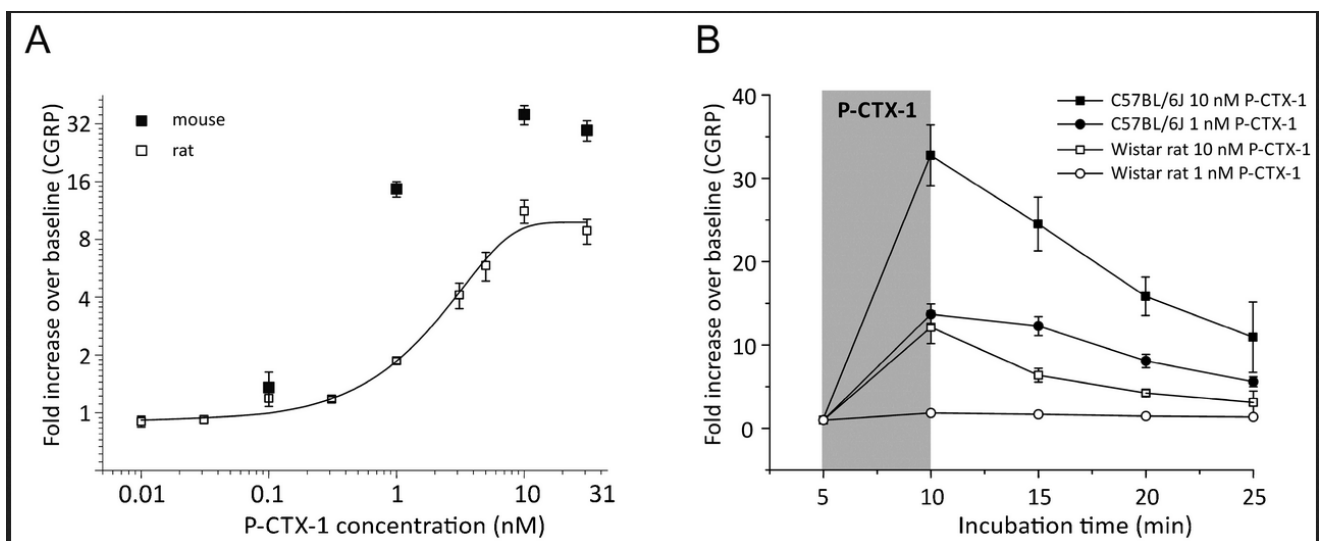


Figure 1: CGRP-release by P-CTX-1 is more effective in mouse than in rat skin. From (2).

- (A) Increasing concentrations of P-CTX-1 augmented the iCGRP release dose-dependently in rat and mouse hindpaw skin. The diagram illustrates fold increase of CGRP as compared to baseline. Error bars represent  $\pm$  SEM;  $n = 4-18$  (see text).
- (B) Release of iCGRP following exposure of rat or mouse hindpaw skin with P-CTX-1 at 1 or 10 nM. Stimulation in mouse tissue was more effective and iCGRP release remained high in the incubation steps after P-CTX-1 was removed. Data are measured as pg/mL and expressed as fold increase over baseline.



# EVALUATION OF CIGUATOXIN-INDUCED CGRP RELEASE IN MOUSE AND RAT SKIN PREPARATIONS

## CONCLUSION

The results of this study show that **P-CTX-1 (Pacific Ciguatoxin-1)** is more effective at inducing **CGRP release in mouse skin as compared to rat skin**, with EC50 concentrations in mouse skin preparations in the low nanomolar range. Since CGRP plays a major role in pain and itch sensing pathways, these findings can help further our understanding of the neurosensorial symptoms of ciguatera, and guide the identification of potential therapeutic targets. **The CGRP Rat ELISA kit from Bertin Bioreagent allows accurate measurements of CGRP in mouse and rat skin preparations.**


## ACKNOWLEDGEMENTS & REFERENCES

Universitätsklinikum  
Erlangen



1. Touska, Filip, et al. "Ciguatoxins evoke potent CGRP release by activation of voltage-gated sodium channel subtypes Nav1. 9, Nav1. 7 and Nav1. 1." *Marine drugs* 15.9 (2017): 269.
2. Zimmermann, K., P. W. Reeh, and B. Averbeck. "S (+)-flurbiprofen but not 5-HT1 agonists suppress basal and stimulated CGRP and PGE2 release from isolated rat dura mater." *Pain* 103.3 (2003): 313-320.
3. Averbeck, B., and P. W. Reeh. "Interactions of inflammatory mediators stimulating release of calcitonin gene-related peptide, substance P and prostaglandin E2 from isolated rat skin." *Neuropharmacology* 40.3 (2001): 416-423.

## DISCOVER OUR COMPREHENSIVE RANGE OF SOLUTIONS

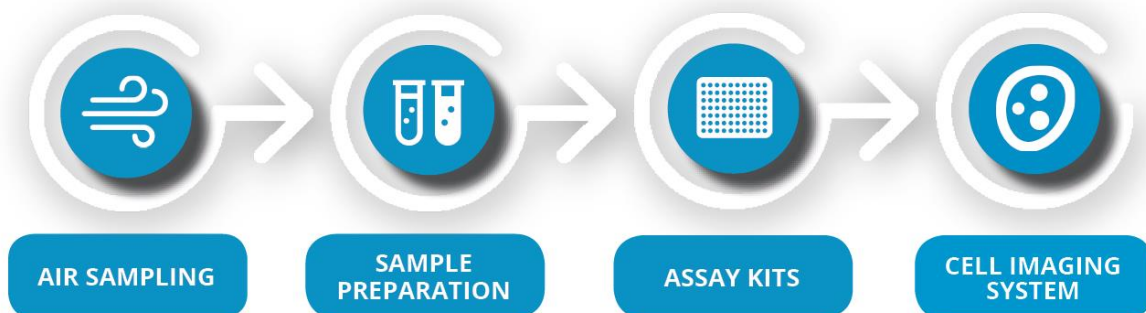

Join the Bertin Instruments community!

1500+ documents available online on the Application center  
[www.bertin-technologies.com/application-center](http://www.bertin-technologies.com/application-center)